

## Environmental Walk-in Chambers



for Stability Testing and Storage





# Environmental Walk-in Chambers for Stability Testing and Storage in Three Different Sizes

WIC	Page 3
Technical Data	Page 4
Tecnical drawings	Page 5
Ordering Options	Page 6

Our Environmental Walk-in Chambers with either 6 m<sup>2</sup>, 9 m<sup>2</sup> or 12 m<sup>2</sup> internal dimensions and 2,3 m height are an excellent choice for stability testing in Pharma, Food & Beverage, Cosmetics, Material Testing and many other fields of application.

Perfect test conditions are maintained in a wide range of temperature from 10°C to 50°C and relative humidity of 20% to 90% RH.

In Pharma, our Walk-in Chambers cover long-term and accelerated stability test conditions as specifed in the ICH Q1A (R2) for new drug substances and drug products.

### **Series WIC**

Our walk-in chambers are designed and produced in Germany, and the temperature and environmental conditioners are made by Parameter Generation and Control in Black Mountain, NC USA. The chambers offer a wide range of features and benefits to perfectly meet requirements in stability testing.



### **PRODUCT FEATURES**

- Relative humidity accuracy of ± 2.5%
- Temperature accuracy of ± 1.5°C
- Duct work for proper distribution of air
- Highly insulated room panels with U value of only 0,21 W/m<sup>2</sup>K
- LED Lighting with motion detector
- High quality room elements, made in Germany
- Interior completely stainless steel
- Exterior service accessibility

### **BENEFITS**


- Stable and uniform environment
- Fully complies with Pharma regulations (ICH)Equipment
- Longevity; corrosion resistant
- Ease of operation
- Environmental friendly refrigerant with GWP value 547
- Maintenance possible with the least possible disruption to the climate conditioning
- Perfectly meets requirements in Pharma industry
- 3 years warranty

Specification:	
Room Elements	Galvanized sheet steel, polyester coated, white, similar to RAL 9010, interior surfaces in stainless steel, 4 inch / 100 mm PU foam U value: $0.21 \text{W/m}^2 \text{K}$
	Floor overlay is chrome-nickel steel, pattern-rolled 5WL
Floor	Surface load: 50.000 N/m². Wheel load (rubber wheel): 1000 N/4 cm²
	Wheel load (steel wheel): Not passable!
	Hinged door (single leaf), door leaf: 4 inch,
	Clear door opening dimensions: 900 mm x 2000 mm
Door	Door version with press lever
	Viewing window, triple insulated, fixed glass
	Clear dimensions of window: 370 mm x 370 mm

### **TECHNICAL DATA**

Description	
Electrical	
Frequency	50 HZ
Voltage (Mains Power)	230 Volts
Phase	1
Full Load Amps	11 Amps, FLA
Typical Running Amps	8 Amps, RLA
Minimum Circuit Ampacity	13 Amps
Maximum External Protector	20 Amps
Maximum Heat of Rejection	5.402 BTU/Hr (1.361 KCAL/Hr)
Heaters, Air	694/918 Watts
Water	
Process Water Max. Pressure	8,6 Bar
Process Water Min. Pressure	0,02 Bar
Process Water Max. Usage	11,4 l/d
Max. Condensate Production	11,4 l/d
Temperature and Relative Environmental Control	
Temperature Range	10-50°C
Acceptable Ambient Range	18-35°C
Temperature Variantion (Spatial)	±1.5 K
Temperature Fluctuation (Temporal)	±0.5 K
Relative Environmental Range	20-93% RH (limited to 7°C dew point)
Relative Humidity Constancy	±2.5 %Rh
Refrigerant	
R450A	0.4 tons C02 / 547 GWP

### TEMPERATURE-ENVIRONMENTAL CHART

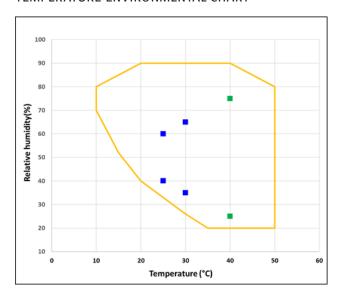

ICH long term stability testing
Accelerated stability testing

### **FEATURES & ACCESSORY**



Lockable door. Even when locked, the door can be opened from the inside at any time.



Lighting with motion sensor.



Interior walls and ceiling: stainless steel
Heavy duty shelves.
Shelf height: 1.800 mm
4 shelf levels standard. Additional shelves optional.
Room lighting with motion detector standard.
Robust stainless steel trolley ramp for safe and convenient access.



Conditioner postitioning at backside of chamber.
Low-noise operation and tidy connection lines.
Technology hidden from front view. External service possible
Mounting of the conditioner on the front is optionally possible.



### **INSTALLATION REQUIREMENTS**

### Installation requirements:


flat, leveled floor electric power supply 230 V / 50 Hz, 11.5 A feed water (tap water or softened water) Drain water hose / connection


in-house transportation to side of installation: by the customer. \\

Chamber: clearance to walls > 40 mm; clearance to wall (door side): 1100 mm


Conditioner: clearance to walls > 915 mm

### WIC 1 | Door on short side

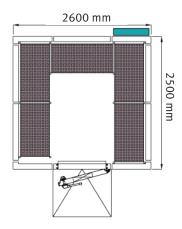


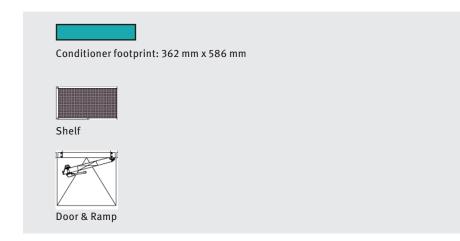


WIC 2 | Door on short side

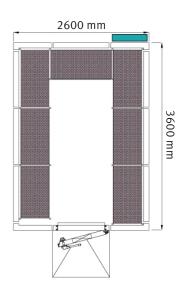


WIC 2 | Door on long side

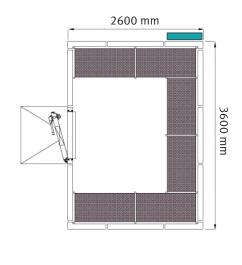

WIC 3 | Door on short side

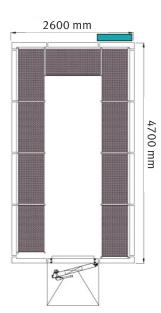

WIC 3 | Door on long side

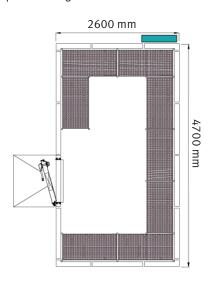

### DIMENSIONS (ROOM W/O CONDITIONER, ACCESSORIES)

Room Sizing Options					
Size	Dimensions	Side 1	Side 2	Height	Usable space (shelves)
WIC 1	outer	2600 mm	2500 mm	2470 mm	187 cb.ft. / 5295 liter
WIC 1	inner	2400 mm	2300 mm	2270 mm	
WIC 2	outer	outer 2600 mm 3600 mm 247	2470 mm	24 -h & /7/74 liber	
WIC 2	inner	2400 mm	3400 mm	2270 mm	— 264 cb.ft. / 7476 liter
WIC 3	outer	2600 mm	4700 mm	2470 mm	2/2 ch ft / 0/05 liter
	inner	2400 mm	4700 mm	2270 mm	— 342 cb.ft. / 9685 liter
*with shelf set					
Conditioner, vertical	outer	362 mm	586 mm	1956 mm	

### ORDERING OPTIONS WALK-IN CHAMBER MODEL

Walk-in Chamber Model	Description	Order No.
WIC 1 (230 V)	Walk-in chamber, size 1	9800-0011
WIC 2 (230 V)	Walk-in chamber, size 2	9800-0012
WIC 3 (230 V)	Walk-in chamber, size 3	9800-0013

### ORDERING OPTIONS ACCESSORY

Accessory	Description	Order No.
Ramp for all WIC models	suitable for all WIC models (230 V)	8012-2239
Shelf set small (for WIC 1 models)	suitable for WIC 1, WIC 1 C (230 V)	8012-2156
Shelf set medium (for WIC 2 models)	suitable for WIC 2, WIC 2 C (230 V)	8012-2157
Shelf set large (for WIC 3 models)	suitable for WIC 3, WIC 3 C (230 V)	8012-2158

### **SERVICES**

Description	Order No.
Installation and Comissioning	DL10-0800
IQ-OQ Qualification	DL40-0200
Temperature & humidity measurement, 13-1 points, at a given climatic point	DL30-0113
Service (Preventive Maintenance)	DL20-1000

## **BINDER** products

CO <sub>2</sub> incubators
Incubators
Cooling incubators
Growth chambers
Ultralow temperature freezers
Drying and heating chambers
Safety drying chambers
Vacuum drying chambers
Constant climate chambers
Dynamic climate chambers
Battery test chambers
Multi Management Software



Never miss a product novelty -

Subscribe and keep up to speed with the BINDER monthly newsletter.

Subscription

go2binder.com/en-newsletter-subscription



Best conditions for your success

BINDER GmbH Im Mittleren Ösch 5 Tuttlingen, Germany

Tel +49 7462 2005 0 info@binder-world.com www.binder-world.com